

Subject Data

ID: ---
Name:
Date of Birth:21-Sep-1677
Age: 340
Gender: ---

Study & Visit Data

Protocol ID: ---
Site ID: ---
Investigator ID: ---
Visit ID: ---
Date & Time: 31-May-2018 14:39:37
Duration: 15:53:22

Machine Data

Device: Savvy 1.0.2 hw=1.23MobECG 1.7.8
Device ID: ---
Software: ---
Manufacturer : ---

Treatment

Referred by:

Analysis Data

Total QRS: 57060
Analyzed Data: 13:48:54, Noise: 01:59:29 (12.5%)

Atrial Fibrillation: -

VPB
Total Number: 149 (0.26%)
Isolated V: 147
Ventricular Couplets: 1
Ventricular Run: 0
V Bigeminy: 0
V Trigeminy: 0
Ventricular Tachycardia: 0

Bradycardia: 0
Pauses: 0

HRV
PNN50: 0.94 %
SDNN: 55.2 ms
RMSSD: 20.1 ms

	Min		Max		Avr
RR (ms)	594	(06:27:41)	1250	(21:59:49)	876
HR (bpm)	59	(22:27:25)	89	(17:15:59)	69

SVPB
Total Number: 21 (0.04%)
Isolated S: 21
Supraventricular Couplets: 0
Supraventricular Run: 0
SV Bigeminy: 0
SV Trigeminy: 0
Supraventricular Tachycardia: 0

Analysis Settings

Analysis Software: CER-S v4.1.0SPRINT1_BETA2 - Beat Detection (ABILE-Pancrazio v1.5.0) - Rhythm Analysis (ABILE-Pancrazio v1.5.0)

BEAT DETECTION

- Leads used for detection: Unknown1

- Noise sensitivity level: High

RHYTHM ANALYSIS

- Prematurity index: 75%

- Prolonged RR Interval Threshold: 150%

- Bradycardia Threshold: 45 bpm, Min. Len.: 10.0 s

- Vent. Tachycardia Thresh.: 150 bpm, Min. Len.: 5.0 s

- Average RR Window size: 8


- Pause Threshold: 2.5 s

- Supravent. Tachycardia Thresh.: 150 bpm, Min. Len.: 5.0 s


Review

Date of Analysis

Reviewer:

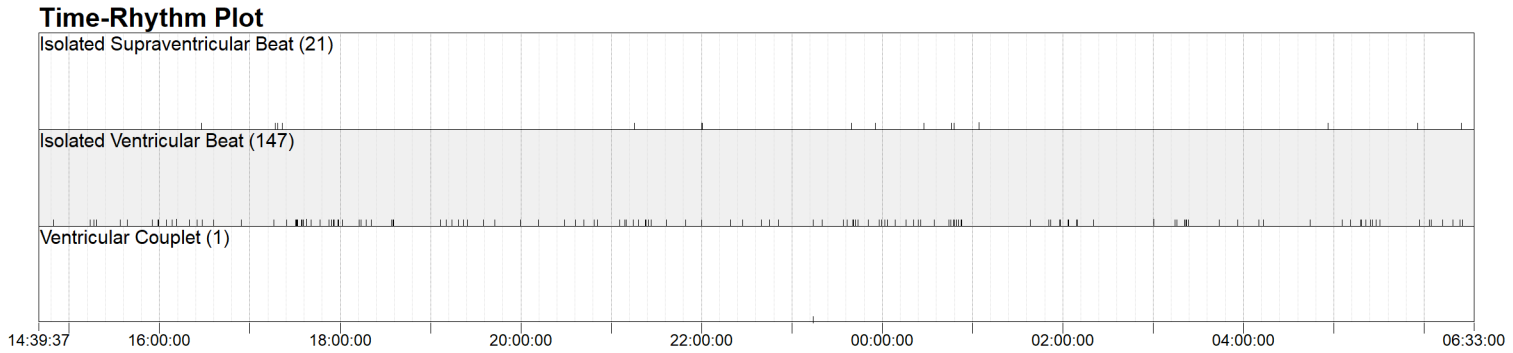

RR Average with extrema and Standard Deviation. (RR window = 10 minute(s))

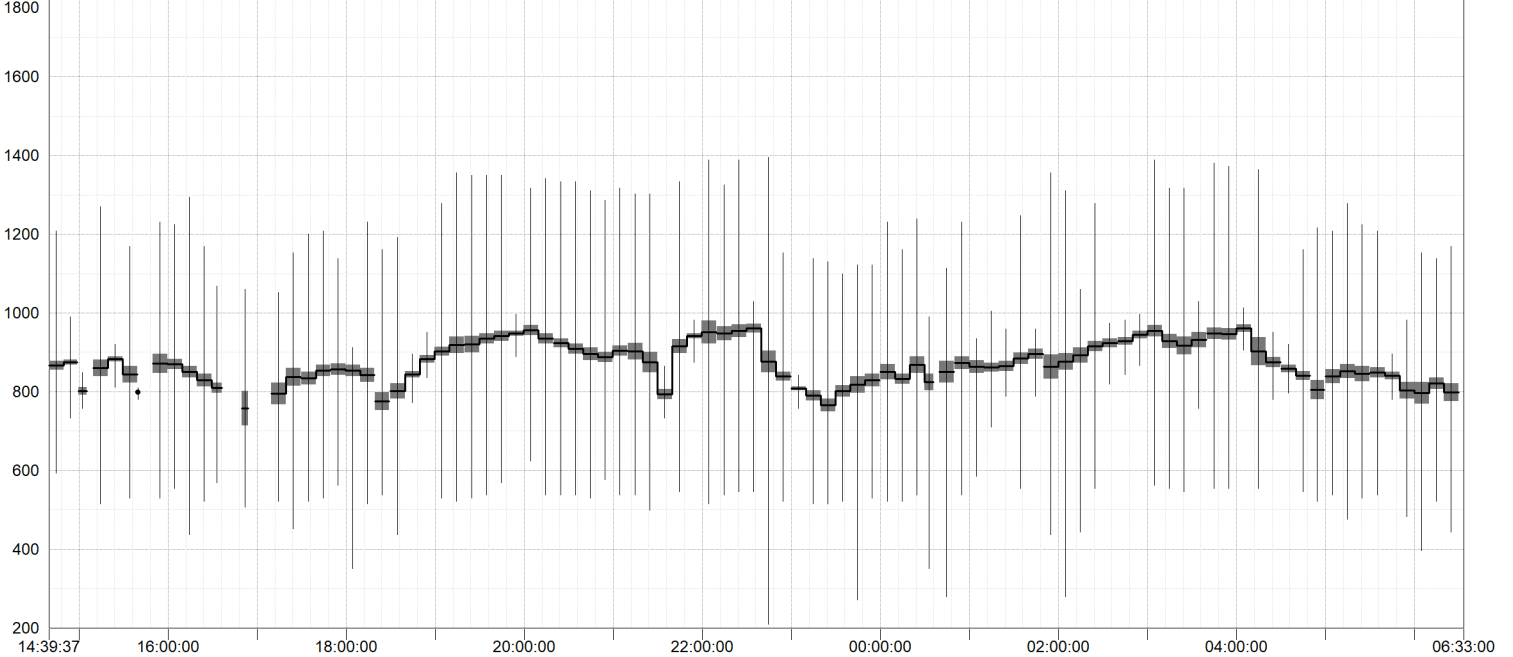

Beats and Rhythm Table - Part 1 of 3

Time	Analyzed Time	RR Min	RR Max	RR Average	RR StDev	HR Min	HR Max	HR Average	HR StDev	Brady	Pauses	Prolonged RR Intervals	SVPB Isolated	SV Couplet	SV Run	SV Bg	SV Tg
14:39:37 - 14:59:59	00:19:25 (95.3%)	734	992	871	26.75	66	72	69	11.04	-	-	-	-	-	-	-	-
15:00:00 - 15:59:59	00:37:01 (61.7%)	734	938	859	38.27	65	80	70	11.13	-	-	-	-	-	-	-	-
16:00:00 - 16:59:59	00:35:26 (59.1%)	742	922	843	44.56	67	83	71	5.79	-	-	-	1	-	-	-	-
17:00:00 - 17:59:59	00:38:39 (64.4%)	633	1125	844	38.73	66	89	71	4.25	-	-	-	3	-	-	-	-
18:00:00 - 18:59:59	00:52:10 (87.0%)	680	992	839	46.42	64	84	72	9.56	-	-	-	-	-	-	-	-
19:00:00 - 19:59:59	00:57:52 (96.5%)	719	1086	930	40.75	61	73	65	8.65	-	-	-	-	-	-	-	-
20:00:00 - 20:59:59	00:59:11 (98.6%)	688	1063	918	38.29	61	70	65	7.29	-	-	-	-	-	-	-	-
21:00:00 - 21:59:59	00:51:16 (85.5%)	734	1250	903	57.04	60	80	67	8.94	-	-	-	4	-	-	-	-
22:00:00 - 22:59:59	00:58:29 (97.5%)	727	1250	921	68.49	59	77	65	8.26	-	-	-	4	-	-	-	-
23:00:00 - 23:59:59	00:52:21 (87.3%)	609	1125	806	30.23	69	87	75	9.09	-	-	-	2	-	-	-	-
00:00:00 - 00:59:59	00:48:35 (81.0%)	680	1164	854	41.17	65	79	70	6.93	-	-	-	3	-	-	-	-
01:00:00 - 01:59:59	00:58:51 (98.1%)	688	1008	874	37.38	63	79	69	10.04	-	-	-	1	-	-	-	-
02:00:00 - 02:59:59	00:58:27 (97.4%)	758	1063	916	53.82	62	75	66	9.02	-	-	-	-	-	-	-	-
03:00:00 - 03:59:59	00:59:53 (99.8%)	758	1031	939	52.27	60	78	64	11.01	-	-	-	-	-	-	-	-
04:00:00 - 04:59:59	00:55:36 (92.7%)	711	1219	875	59.98	60	83	69	10.70	-	-	-	1	-	-	-	-
05:00:00 - 05:59:59	00:57:11 (95.3%)	641	922	840	43.18	67	84	71	8.60	-	-	-	1	-	-	-	-
06:00:00 - 06:32:59	00:28:23 (86.1%)	594	898	808	54.97	69	88	74	7.22	-	-	-	1	-	-	-	-
All Days: 08:00:00 - 20:00:00	04:00:36 (75.1%)	633	1125	867	56.77	61	89	69	7.93	-	-	-	4	-	-	-	-
All Nights: 23:00:00 - 06:00:00	06:30:56 (93.1%)	609	1219	872	67.98	60	87	69	4.55	-	-	-	8	-	-	-	-
Entire Recording	13:48:54 (86.9%)	594	1250	876	59.43	59	89	69	10.42	-	-	-	21	-	-	-	-

Beats and Rhythm Table - Part 2 of 3

Time	Analyzed Time	SV Tachy	ATachy	AFib	AFlutter	V Isolated	V Couplet	V Run	V Bg	V Tg	V Tachy	VFib	VFlutter	Torsade de Pointes	AVB	N Beats	V Beats
14:39:37 - 14:59:59	00:19:25 (95.3%)	-	-	-	-	1	-	-	-	-	-	-	-	-	-	1338	1
15:00:00 - 15:59:59	00:37:01 (61.7%)	-	-	-	-	8	-	-	-	-	-	-	-	-	-	2599	8
16:00:00 - 16:59:59	00:35:26 (59.1%)	-	-	-	-	8	-	-	-	-	-	-	-	-	-	2527	8
17:00:00 - 17:59:59	00:38:39 (64.4%)	-	-	-	-	20	-	-	-	-	-	-	-	-	-	2751	20
18:00:00 - 18:59:59	00:52:10 (87.0%)	-	-	-	-	10	-	-	-	-	-	-	-	-	-	3753	10
19:00:00 - 19:59:59	00:57:52 (96.5%)	-	-	-	-	9	-	-	-	-	-	-	-	-	-	3739	9
20:00:00 - 20:59:59	00:59:11 (98.6%)	-	-	-	-	6	-	-	-	-	-	-	-	-	-	3870	6
21:00:00 - 21:59:59	00:51:16 (85.5%)	-	-	-	-	12	-	-	-	-	-	-	-	-	-	3410	12
22:00:00 - 22:59:59	00:58:29 (97.5%)	-	-	-	-	5	-	-	-	-	-	-	-	-	-	3815	5
23:00:00 - 23:59:59	00:52:21 (87.3%)	-	-	-	-	13	1	-	-	-	-	-	-	-	-	3912	15
00:00:00 - 00:59:59	00:48:35 (81.0%)	-	-	-	-	16	-	-	-	-	-	-	-	-	-	3446	16
01:00:00 - 01:59:59	00:58:51 (98.1%)	-	-	-	-	5	-	-	-	-	-	-	-	-	-	4042	5
02:00:00 - 02:59:59	00:58:27 (97.4%)	-	-	-	-	5	-	-	-	-	-	-	-	-	-	3836	5
03:00:00 - 03:59:59	00:59:53 (99.8%)	-	-	-	-	9	-	-	-	-	-	-	-	-	-	3819	9
04:00:00 - 04:59:59	00:55:36 (92.7%)	-	-	-	-	3	-	-	-	-	-	-	-	-	-	3827	3
05:00:00 - 05:59:59	00:57:11 (95.3%)	-	-	-	-	10	-	-	-	-	-	-	-	-	-	4091	10
06:00:00 - 06:32:59	00:28:23 (86.1%)	-	-	-	-	7	-	-	-	-	-	-	-	-	-	2115	7
All Days: 08:00:00 - 20:00:00	04:00:36 (75.1%)	-	-	-	-	56	-	-	-	-	-	-	-	-	-	16707	56
All Nights: 23:00:00 - 06:00:00	06:30:56 (93.1%)	-	-	-	-	61	1	-	-	-	-	-	-	-	-	26973	63
Entire Recording	13:48:54 (86.9%)	-	-	-	-	147	1	-	-	-	-	-	-	-	-	56890	149


Beats and Rhythm Table - Part 3 of 3


Time	Analyzed Time	S Beats
14:39:37 - 14:59:59	00:19:25 (95.3%)	-
15:00:00 - 15:59:59	00:37:01 (61.7%)	-
16:00:00 - 16:59:59	00:35:26 (59.1%)	1
17:00:00 - 17:59:59	00:38:39 (64.4%)	3
18:00:00 - 18:59:59	00:52:10 (87.0%)	-
19:00:00 - 19:59:59	00:57:52 (96.5%)	-
20:00:00 - 20:59:59	00:59:11 (98.6%)	-
21:00:00 - 21:59:59	00:51:16 (85.5%)	4
22:00:00 - 22:59:59	00:58:29 (97.5%)	4
23:00:00 - 23:59:59	00:52:21 (87.3%)	2
00:00:00 - 00:59:59	00:48:35 (81.0%)	3
01:00:00 - 01:59:59	00:58:51 (98.1%)	1
02:00:00 - 02:59:59	00:58:27 (97.4%)	-
03:00:00 - 03:59:59	00:59:53 (99.8%)	-
04:00:00 - 04:59:59	00:55:36 (92.7%)	1
05:00:00 - 05:59:59	00:57:11 (95.3%)	1
06:00:00 - 06:32:59	00:28:23 (86.1%)	1
All Days: 08:00:00 - 20:00:00	04:00:36 (75.1%)	4
All Nights: 23:00:00 - 06:00:00	06:30:56 (93.1%)	8
Entire Recording	13:48:54 (86.9%)	21

56890 Beats in 84 Templates [page 1 of 4]

Templates N


50 mm/s 10 mm/mV

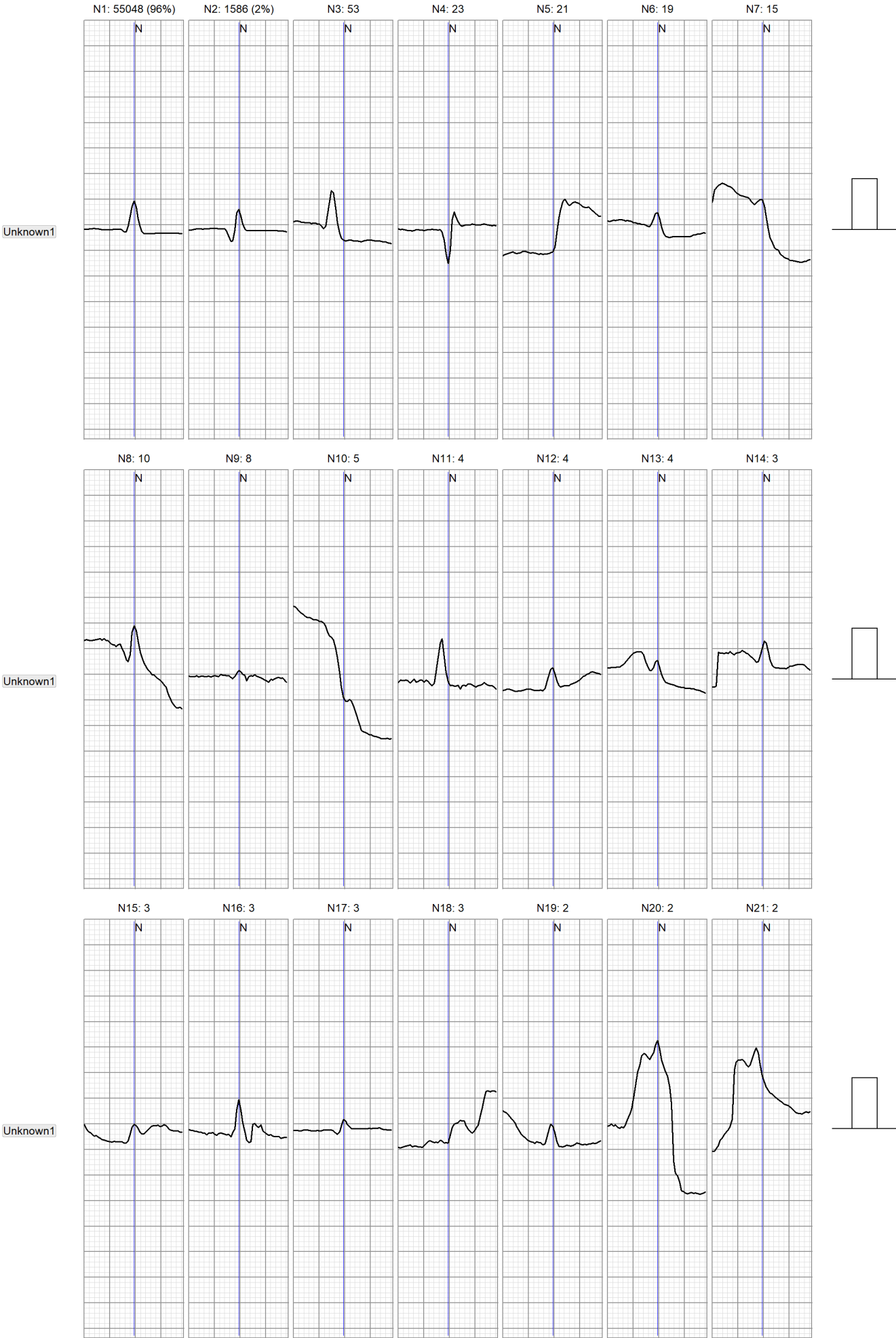

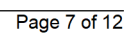


Templates N

50 mm/s 10 mm/mV

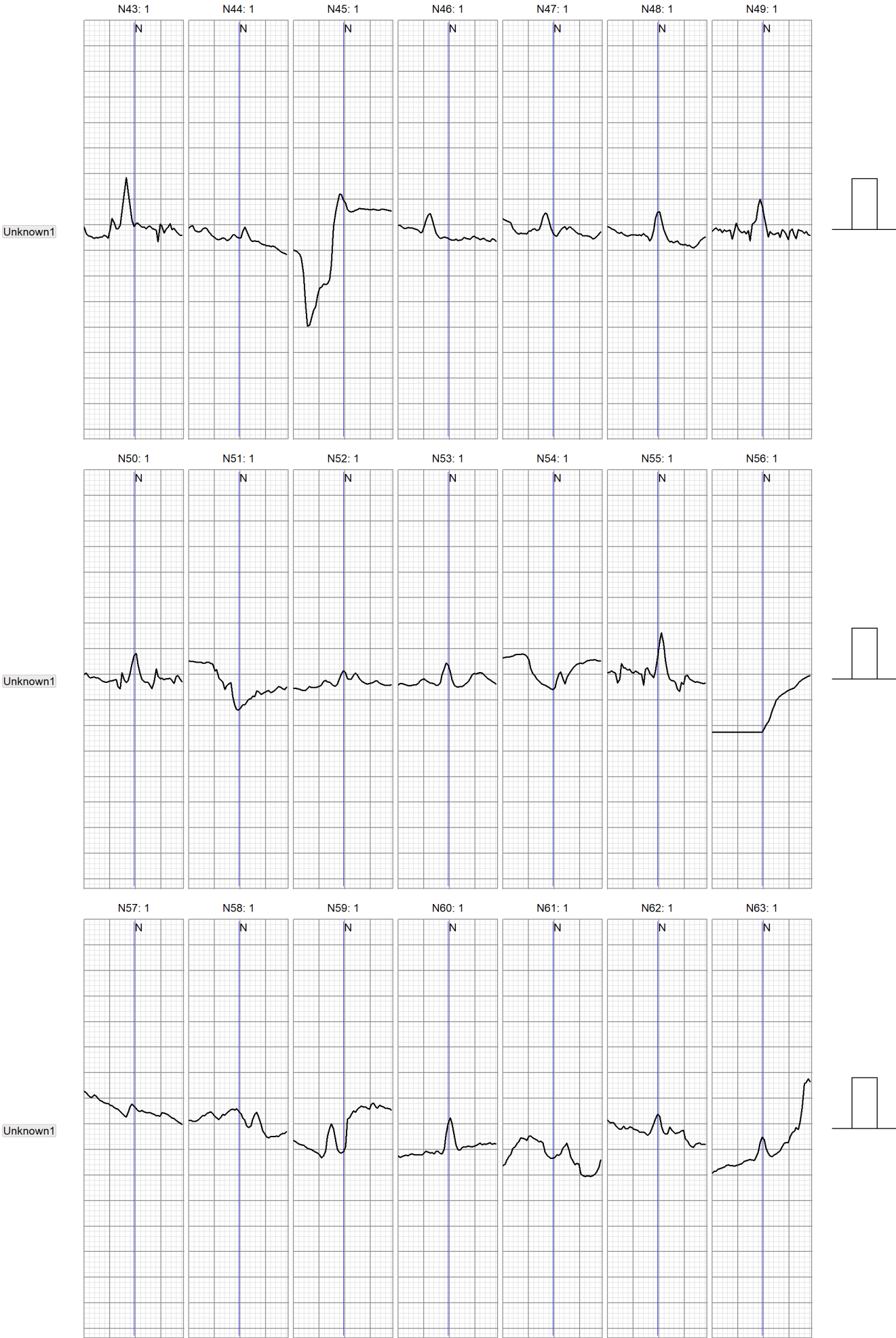


Templates N

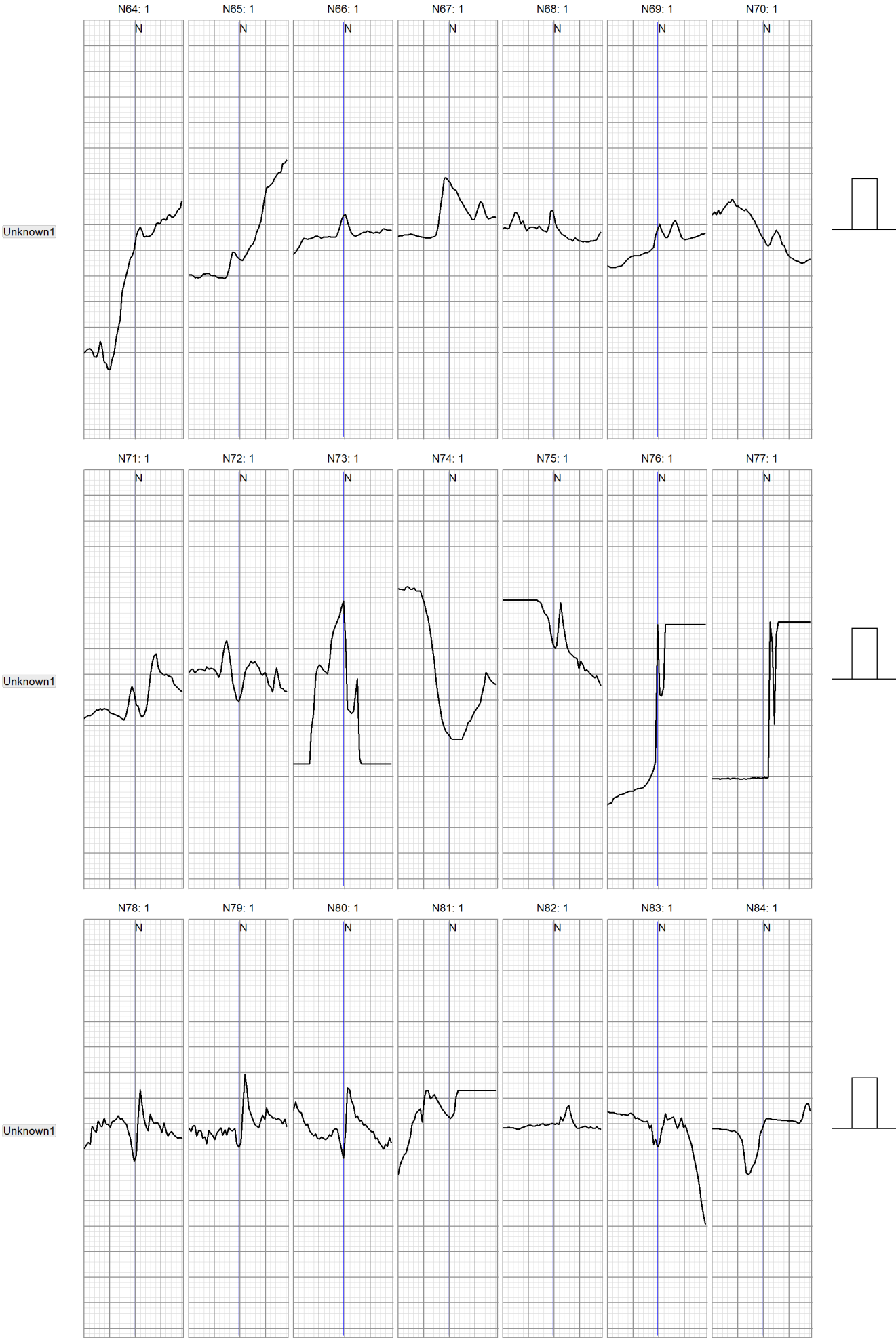

149 Beats in 15 Templates [page 1 of 1]

Templates V

50 mm/s 10 mm/mV

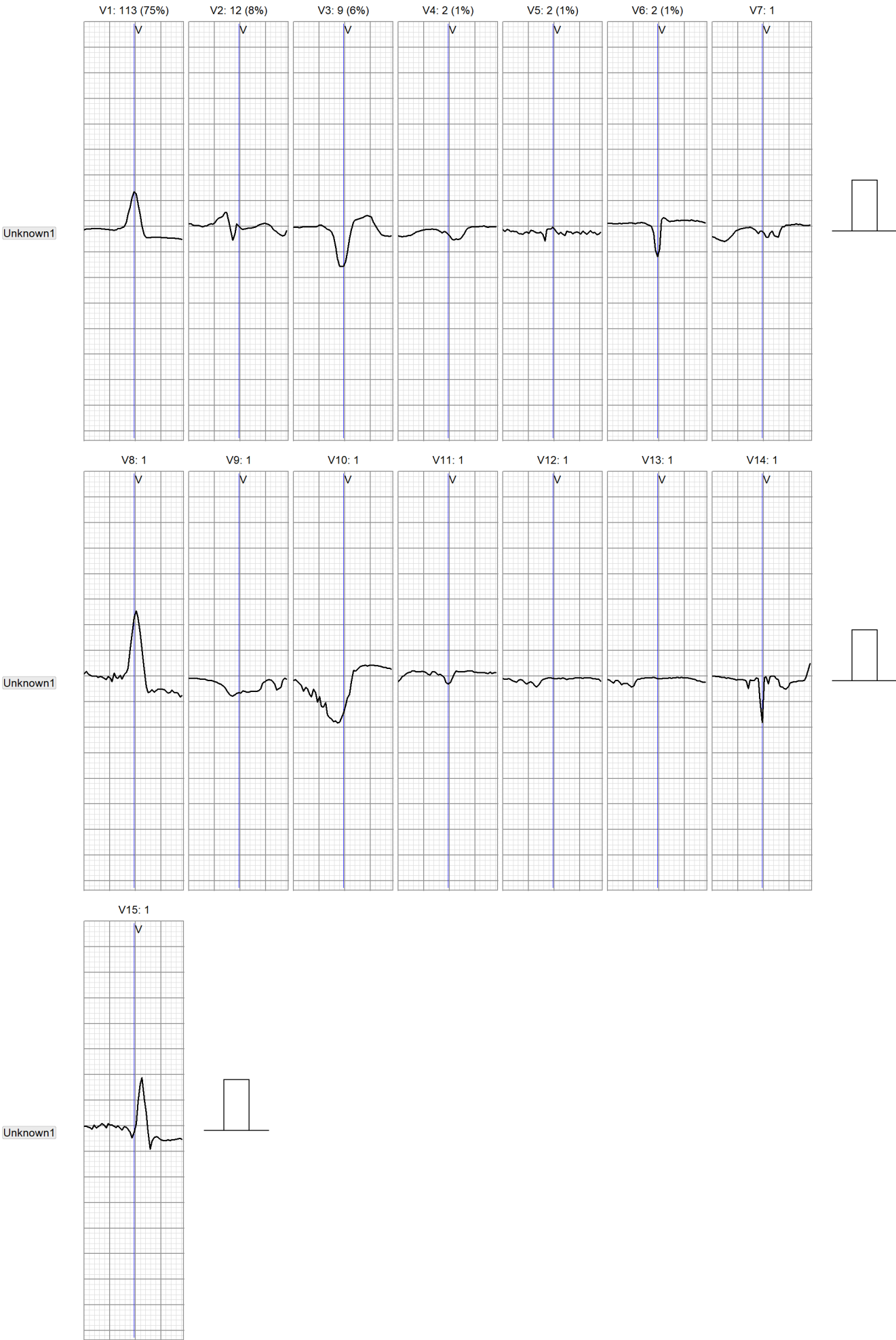


21 Beats in 6 Templates [page 1 of 1]

Templates S

50 mm/s 10 mm/mV

Unknown1

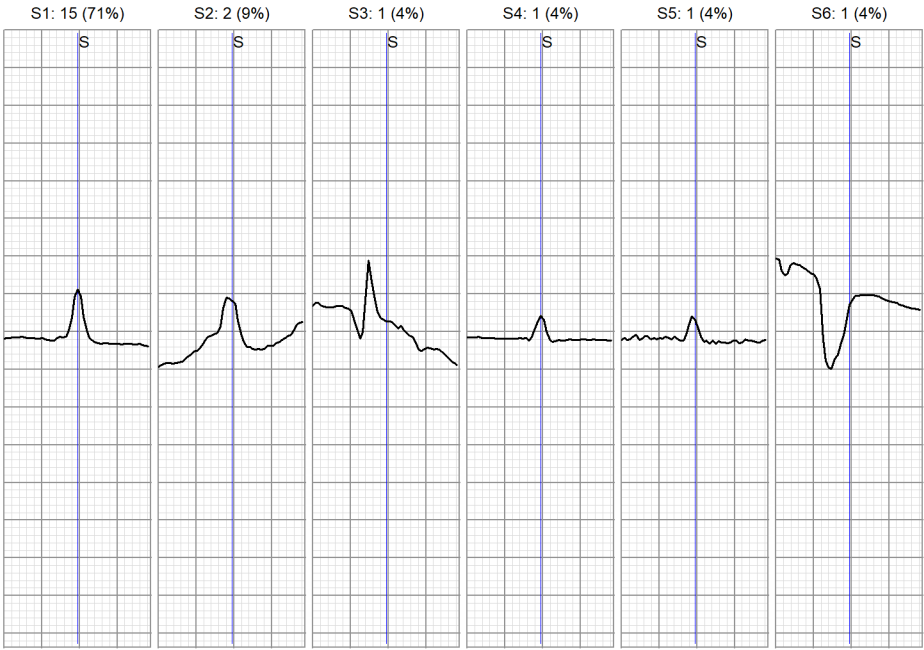


Automatic Strips

Isolated Supraventricular Beat (Shortest): starting at 01-Jun-2018 06:24:09 (7 sec)


25 mm/s 5 mm/mV

**Isolated Ventricular Beat (Shortest):** starting at 31-May-2018 15:13:48 (7 sec)

25 mm/s 5 mm/mV

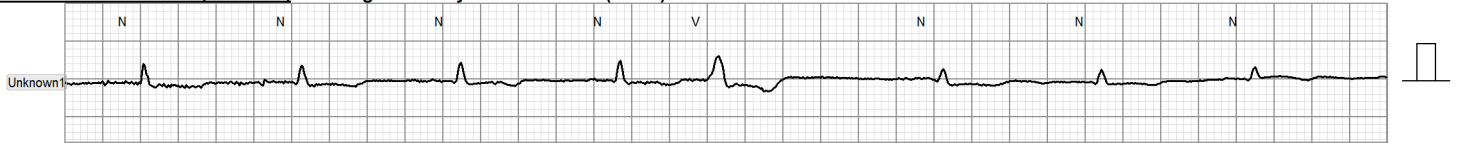**Ventricular Couplet (Shortest):** starting at 31-May-2018 23:13:58 (7 sec)


25 mm/s 5 mm/mV

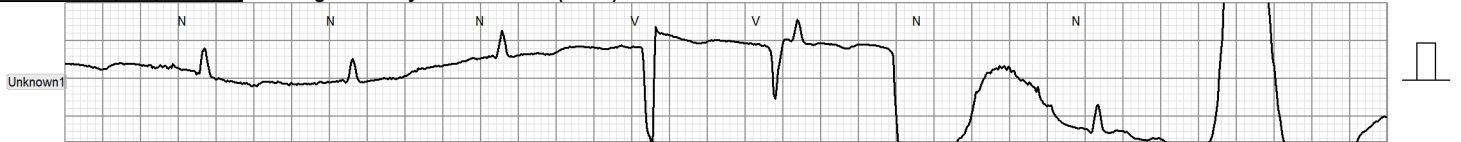**Ventricular Beat of template V1:** starting at 31-May-2018 14:49:25 (7 sec)


25 mm/s 5 mm/mV

**Ventricular Beat of template V2:** starting at 31-May-2018 17:30:13 (7 sec)


25 mm/s 5 mm/mV

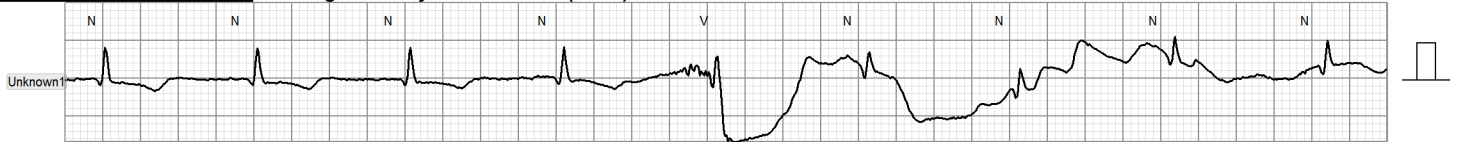**Ventricular Beat of template V3:** starting at 31-May-2018 18:34:49 (7 sec)


25 mm/s 5 mm/mV

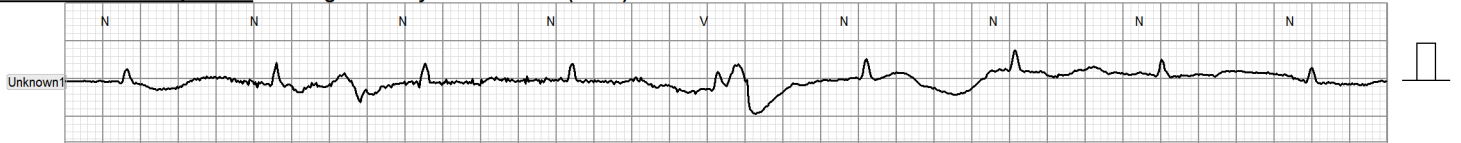**Ventricular Beat of template V4:** starting at 31-May-2018 16:11:24 (7 sec)

25 mm/s 5 mm/mV

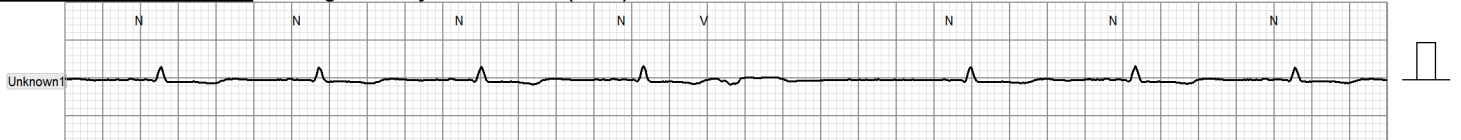**Ventricular Beat of template V5:** starting at 31-May-2018 22:45:04 (7 sec)

25 mm/s 5 mm/mV

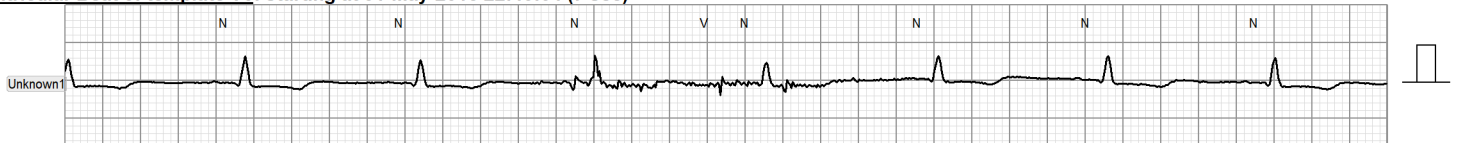**Ventricular Beat of template V6:** starting at 31-May-2018 17:15:51 (7 sec)

25 mm/s 5 mm/mV

